



Yes, you can.[®]

Invacare[®] Pronto M61

A fresh approach to mobility

The sleek and stylish *Invacare[®] Pronto M61* is an affordable centre wheel drive power chair that features *Invacare[®]'s Sure Step[®] Technology* for smoother driving over lower thresholds, and with 6 wheels on the ground it ensures optimum safety & stability. Designed around a compact base, the **Pronto M61** provides excellent manoeuvrability especially in confined spaces, and the electric seat provides 12 cm of elevation.



Optimum Comfort

The padded, shaped van seat ensures you are comfortable for day long use, while the popular configurable seating system offers a wide range of seat sizes allowing you to achieve an optimum fit which aids postural management. The simple to operate seat will electrically raise you up to 12 cm. In the elevated position the **Invacare® Pronto M61** provides increased accessibility, which is ideal in an office environment or when you are out and about. For easier transfers into or out of the **Pronto M61**, the seat will rotate 90 degrees in either direction, and the van seat arm rest can flip up or alternatively the armrests can be removed altogether.



Upwardly mobile

With its centre wheel drive and its compact dimensions, the **Pronto M61** gives excellent manoeuvring in very small spaces and has a very tight turning circle. **Invacare®'s Sure Step®** technology gives a smooth ride over low thresholds of up to 5 cm while maintaining optimum stability at all times. The powerful motors lift the front castors to assist changes in height for a smoother, less jarring ride over rough surfaces.



Tension adjustable upholstery



Comfort seat




Comfort back



Padded seat

Technical data

Van Seat 45,5 - 50 cm	41 - 46 cm	47 - 59 cm*1	0 - 30°	61 cm	68 cm ³	89 - 94 cm
Config Seat 33 - 43 cm	32 - 50 cm	50 - 62 cm ²	0 - 12°	(excluding joystick)	97 - 102 cm ⁴	96 - 110 cm
R = 48 cm	5 cm	9° (16%)	2 x 30 Ah	2 x 180 W	16 km	6,4 km/h
		*1 Seat height is measured including cushion, with 4° tilt angle *2 Seat height is measured including cushion, with 9° tilt angle *3 with footboard up *4 with footboard down, with side mounted legrests *5 Drive wheels (flat free) *6 Front and rear casters				
70 kg (91 kg with batteries)	136 kg					

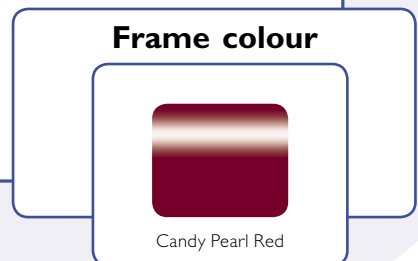

All information quoted is believed to be correct at time of print. Invacare® reserves the right to alter product specifications without prior consultation.



INVACARE® Limited
 South Road · Bridgend Industrial Estate · BRIDGEND CF31 3PY
 Tel. +44 (0) 1656 647327 · Fax +44 (0) 1656 649016
 E-mail: uk@invacare.com · www.invacare.co.uk
 Sales Order E-mail: ordersuk@invacare.com

INVACARE® Ireland Ltd
 Unit 5 · Seatown Business Campus
 Seatown Road · SWORDS · County Dublin
 Tel. (+353) 181 107 84 · Fax (+353) 181 107 85
 Email: ireland@invacare.com · www.invacare.ie

Pronto M61 - UK - 09/07 - Part-No. 1468372
 © Invacare® International SARL 2007. All rights reserved.

