

Uniclass L72101:N5352	EPIC H122:Y311
CI/SfB (74.22)	X (U35)

EB/05/2004

E A S I B A T H

Assisted bathing, showering and changing



Kingkraft
Kingkraft

EASIBATH



Assisted bathing, showering and changing



Easibath was introduced to address the problem of high dependency bathing. It acts as a variable height bath, reducing the risk of back strain, while also serving as a showering platform and changing bench. The load bearing fold side has the twin benefits of allowing ambulant users to enter the bath, as well as giving a Carer access for changing and dressing the user. For non-ambulant users, Easibath is designed for use with fixed, mobile or ceiling mounted hoist systems.

Easibath is available in a range of standard sizes and can also be specially manufactured to the customers' own specification, including Special Order units which are both longer and wider than standard models. A base and side cushion, along with a choice of backrest (not required with 1700 model) are supplied with the bath as standard and a comprehensive range of supports and mouldable vacuum cushions can be used to customise the support provided. Please see our cushions and supports literature for full details.

A combination of supports for the needs of the individual bather is normally determined by demonstration of the bath in the user's home or other convenient location. Kingkraft understand that, over time, the needs of the bather may change and we are available to give advice and supply alternative supports to ensure that the most suitable position for the bather is maintained.

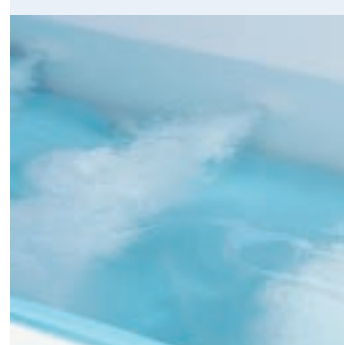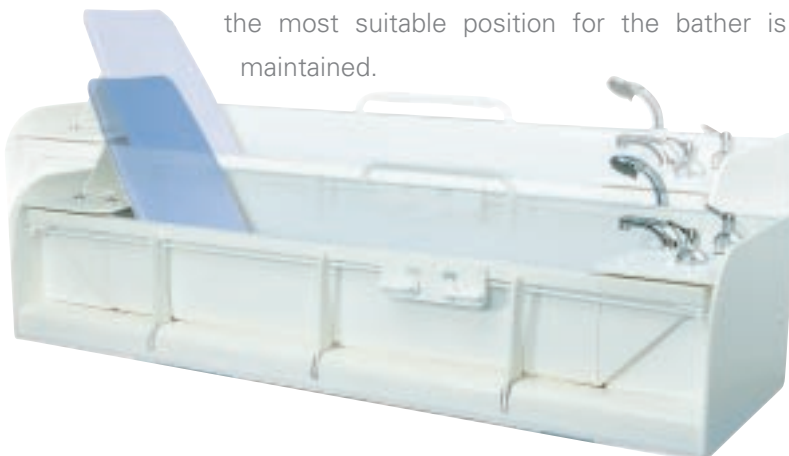


standard features

- Load bearing dropside design
- Raise and lower
- Side and base cushion
- Choice of backrest with cushion
- Integral shower unit
- Fixed temperature thermostatic mixer valves
- Chrome taps and spout
- Left or right handed tap position

additional options

- Whirlpool system
- Wall mounted adjustable mixer valves
- Twin waste
- Custom support systems
- Custom size manufacture
- Fixed height model
- Grab handles

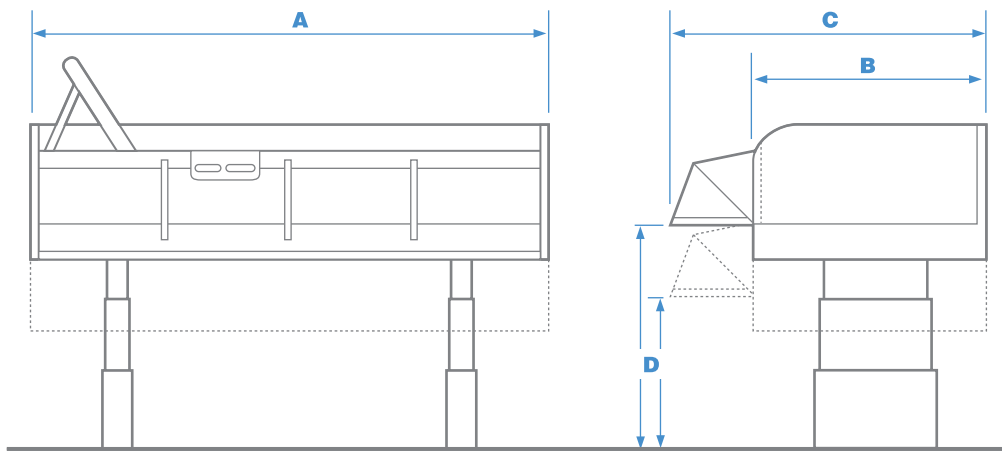

Model No.	A Length	B Width side up	C Width side down	D Base height approx	Internal length	Internal width
VARIABLE HEIGHT MODEL						
1700	1700	885	1065	460-860 variable	1285	645
1900	1837	760	910		1292	525
2100	2116	760	940		1525	525
2100EW	2116	910	1090		1525	675
FIXED HEIGHT MODEL						
1700	1700	885	1065	Fixed between: Min. 200 Max. 860	1285	645
1900	1837	760	910		1292	525
2100	2116	760	940		1525	525
2100EW	2116	910	1090		1525	675

All dimension expressed in millimetres

Maximum Load : 130kg/21st

Kingkraft Customer Care

The following services are offered to make the process of choosing your Kingkraft product as straight forward as possible. Kingkraft Care continues following delivery to ensure the products continue to provide the maximum benefit.

- Product demonstration and assessments
- Bathroom layout design and drawing
- Technical support both prior to order and following delivery
- Product commissioning
- 12 month product warranty (extended warranties available)
- On-going product servicing and preventative maintenance schemes

For full details of these services please contact our customer service department or refer to our product information folder.



Kingkraft Limited, 26D Orgreave Crescent, Dore House Industrial Estate, Sheffield S13 9NQ
 Telephone: 0114 269 0697 Fax: 0114 269 5145 Email: info@kingkraft.co.uk Web: www.kingkraft.co.uk

Due to the process of continuous development we reserve the right to alter specification without notice.