

Colt TWIN

HIGH-PERFORMANCE OPERATION
AND ALL-NEW FEATURES AT A GREATER
VALUE THAN EVER BEFORE

COLT TWIN

The sleek, sporty **Colt Twin** scooter delivers high performance operation and features unique dual front wheels to deliver the maximum stability of a four wheel scooter with the superb manoeuvrability of a three wheel scooter. With a full complement of advanced standard features like exclusive low-profile tyres, wraparound delta tiller, and feather-touch the Colt Twin is a greater value than ever before.



BUILT WITH PRIDE

Pride
Mobility Products Ltd.




(Tel) 01869 324 600 • (Fax) 01869 323 070 • www.pridemobility.com

Colt Twin

Weight capacity:	25 Stone (159 kg)
Maximum speed ¹ :	Up to 6 km/h (4 mph)
Ground clearance ² :	8.25 cm (3.25")
Turning radius:	107 cm (42")
Overall length ² :	117 cm (46") at wheels
Overall width ² :	56.5 cm (22.25")
Range per charge ¹ :	Up to 24 km (15 miles)
Tyre Size:	Front: Dual 9 cm x 23 cm (3.5" x 9") solid Rear: 9 cm x 25 cm (3.5" x 10") solid
Total Weight Without Batteries:	58 kg (128.5 lbs.)
Heaviest piece when disassembled:	Front: 25 kg (55 lbs.)
Standard seating:	Type: Foldable molded plastic Weight: 12.5 kg (27.5 lbs.) Material: Black Vinyl Dimensions: 46 cm (18") width 43 cm (17") depth
Seat-to-Deck Range:	44.5 cm - 48.25 cm (17.5" - 19")
Drive System:	Rear-wheel drive, sealed transaxle, 24-volt DC motor
Dual Braking System:	Electronic, regenerative, and electromechanical
Battery Requirements ³ :	Two 12-volt, deep cycle Size: (2) 30 amp Weight: 14.5 kg (32.5 lbs.) each
Battery Charger:	5-amp, on-board

Colours:

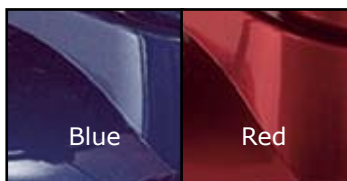

¹Varies with user weight, terrain type, battery amp-hour (AH), battery charge, battery condition and tyre condition.

²Due to manufacturing tolerances and continual product improvement, this specification can be subject to a variance of (+ or - 3%).

³AGM or Gel-Cell type required.

NOTE: All specifications subject to change without notice.

Features

- Dual front wheels for added stability
- Exclusive one-hand Feather-touch disassembly provides the ultimate quick transport system
- Auto-connecting front-to-rear harness
- Delta tiller with wraparound handles
- Pride's exclusive stylish, light weight, non-marking, low profile wheels
- Backlit battery gauge
- Long-lasting, high-intensity LED headlight offers optimal pathway illumination
- Lightweight, seat with viscoelastic foam inserts for increased comfort
- Includes pinchless hinges & eight-position seat rotation
- Easy-grip tiller adjustment knob
- Easy access tiller-mounted charger port for standard dual voltage off-board charger

Options

- Forearm crutch holder
- Cane/crutch holder
- Oxygen tank holder
- Cup holder
- Rear basket
- Safety flag
- Walker holder
- Weather cover
- Saddlebag



The information contained herein is correct at the time of publication; due to the manufacturers commitment to constant improvement and development, they reserve the right to alter specifications without prior notice. Theoretical range is calculated under test conditions in accordance with European standards. Actual range in normal use will depend on many factors, including the condition of the vehicle and its batteries, the weight of the driver, correct tyre pressure, ambient temperature, and the gradient and surface of the road or pavement.

32 Wedgwood Road, Bicester, Oxfordshire OX26 4UL
(Tel) 01869 324 600 • (Fax) 01869 323 070 • www.pridemobility.com

Pride
Mobility Products Ltd.